

Characterization of the small RNA transcriptome using the bebionextgen python framework

Lorena Pantano @lopantano Harvard TH Chan School of Public Health

SMAIL KNA RNA molecules of 18-36 nts long with regulation function

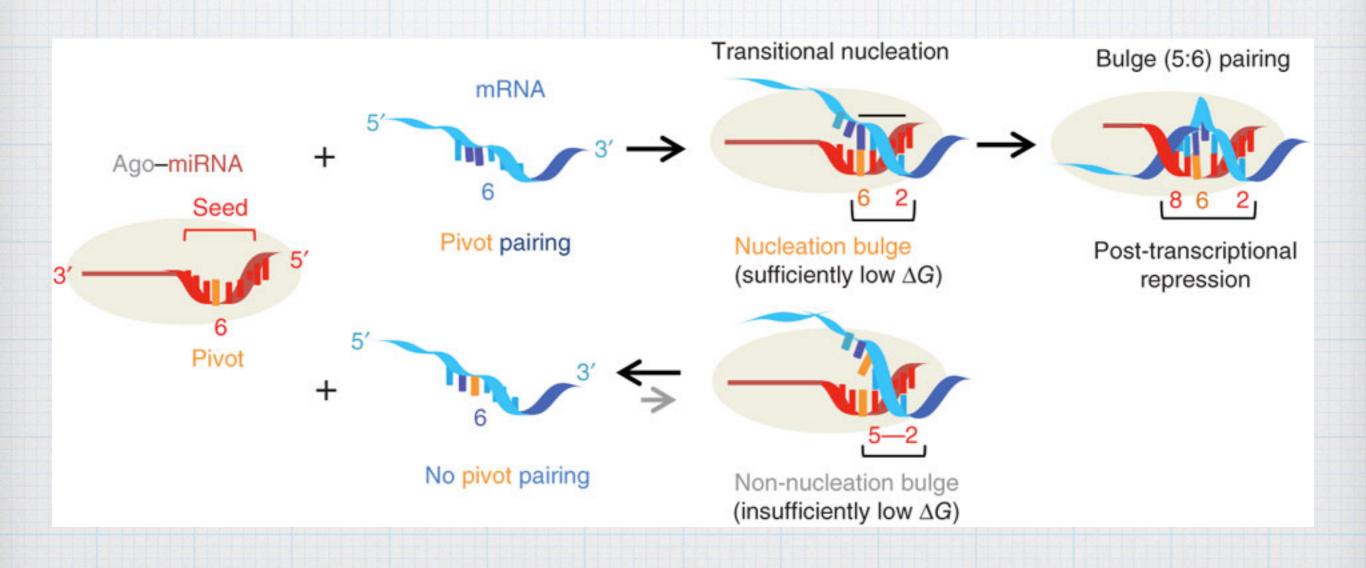

MiKNA axon outgrowth

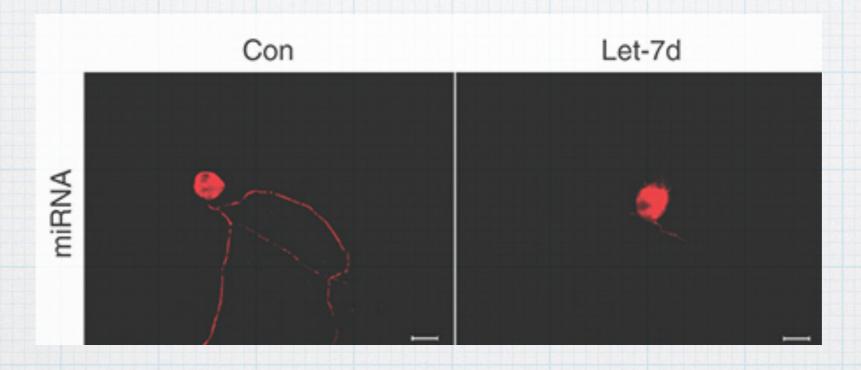


Let-7 microRNAs Regenerate Peripheral Nerve Regeneration by Targeting Nerve Growth Factor Shiying Li, Xinghui Wang, Yun Gu, Chu Chen, Yaxian Wang, Jie Liu, Wen Hu, Bin Yu, Yongjun Wang, Fei Ding, Yan Liu and Xiaosong Gu

isomiks


hsa-miR-24-1-5pGGUGCCUACUGAGCUGAUAUC	hsa-miR-24-3p
GUGCCUACUGAGCUGAUAUCAGU	
GUGCCUACUGAGCUGAUAUCAG	
<u>UGCCUACUGAGCUGAUAUCA</u>	
<u>UGCCUACUGAGCUGAUAUC</u> <u>UGCCUACUGAGCUGAUA</u>	
CCUACUGAGCUGAUAUCA	

types of isomiks

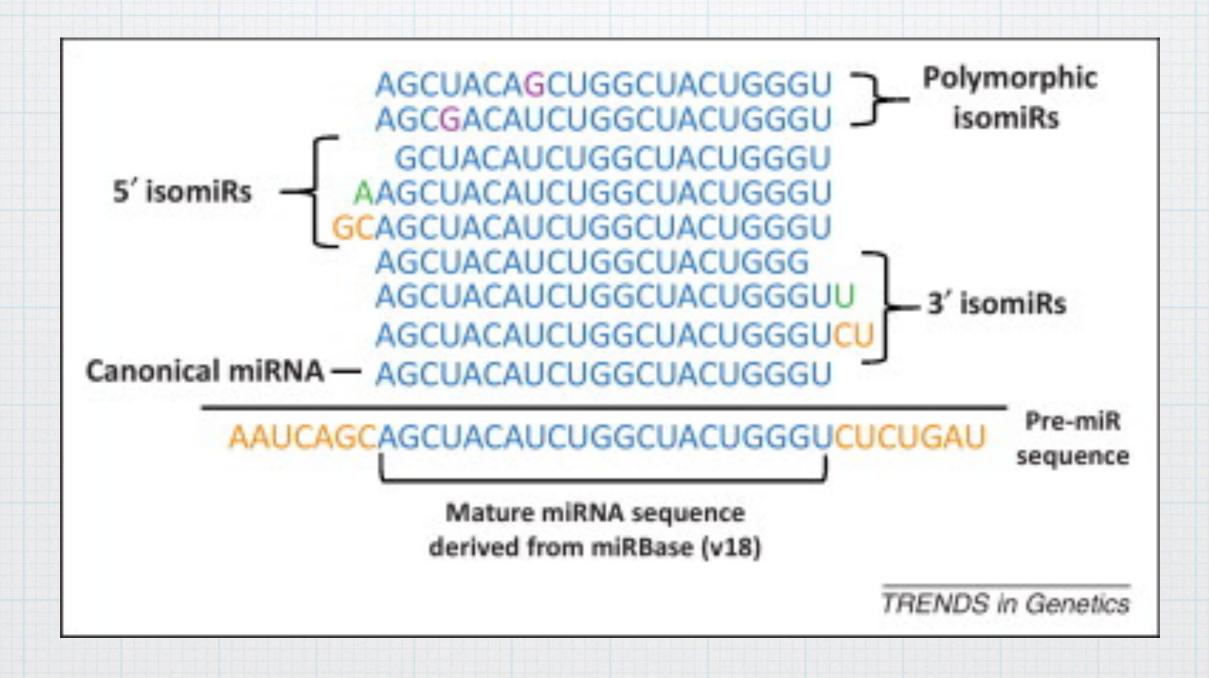


isomiks

Gene expression

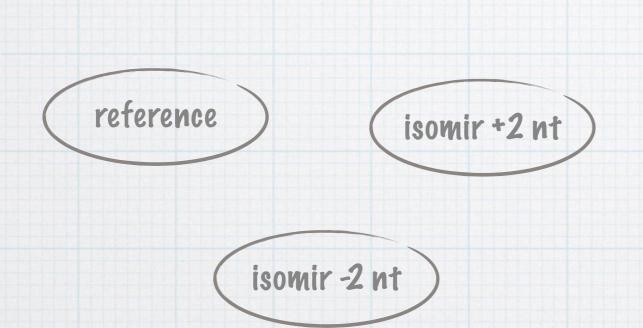


transfected mammary cells line derived from metastatic site

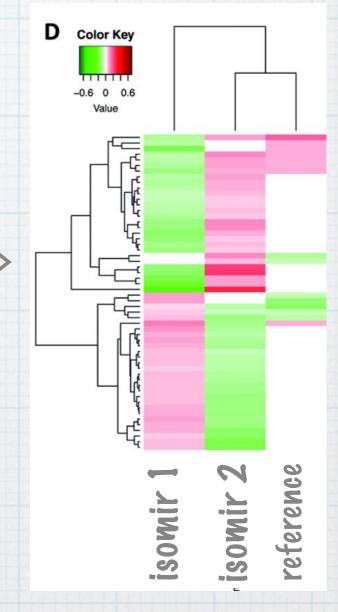


Aristeidis G. Telonis et al. Nucl. Acids Res. 2015;nar.gkv922

small tRNAs

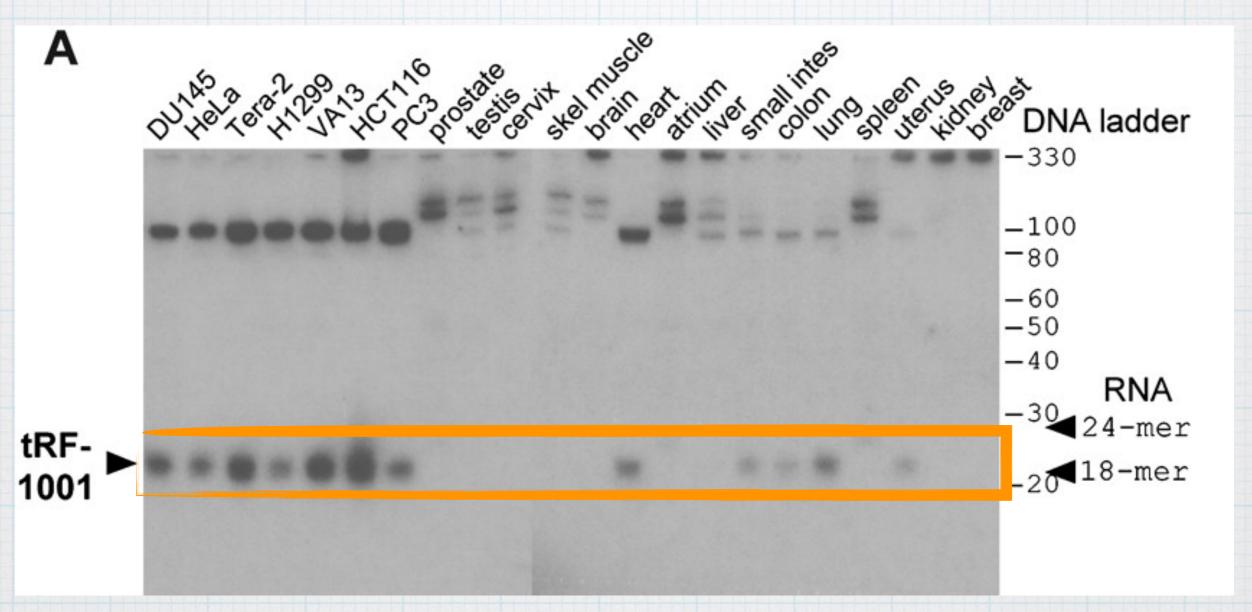

Yong Sun Lee et al. Genes Dev. 2009;23:2639-2649

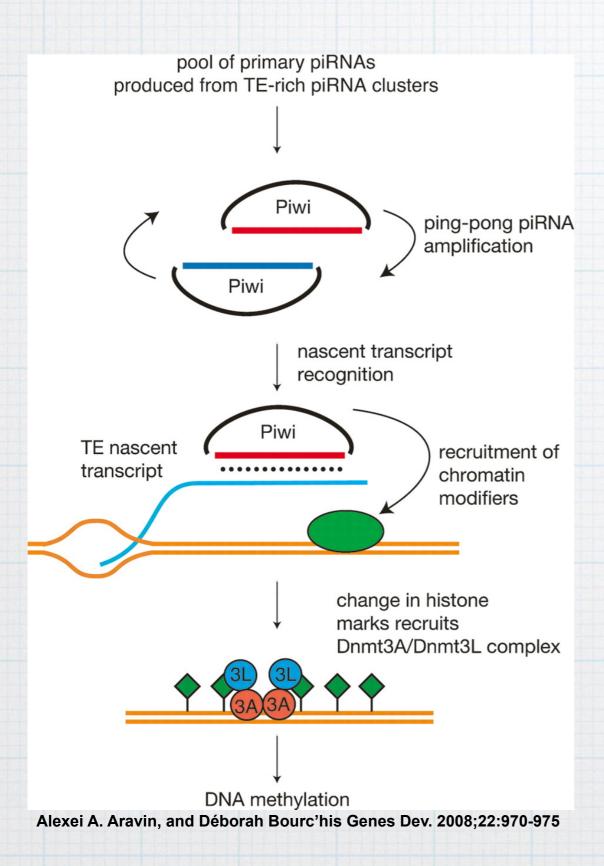


piRNAs

challenges


- * isomiRs detection
- * small RNAs coming from multiple precursors over the genome (multi-mapped reads can be 40% of the data.)
- * differentiate degradation and functional molecules
- * non-model organism

bcbio-nextgen

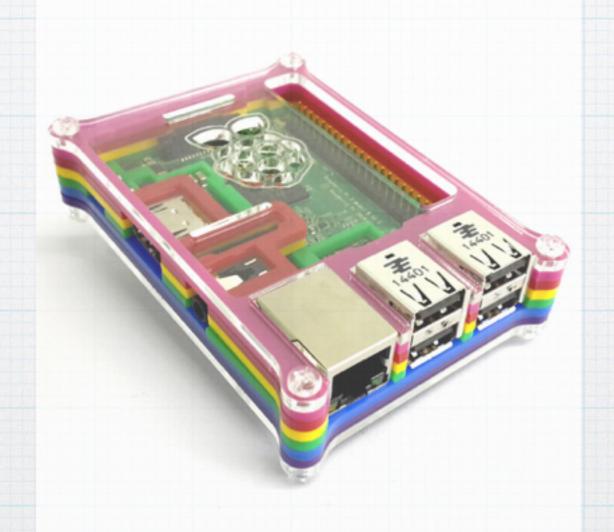

Variant calling, RNA-seq, small RNA-seq over 200 peer reviewed tools BIOCONDA*

Period: 1 month ▼

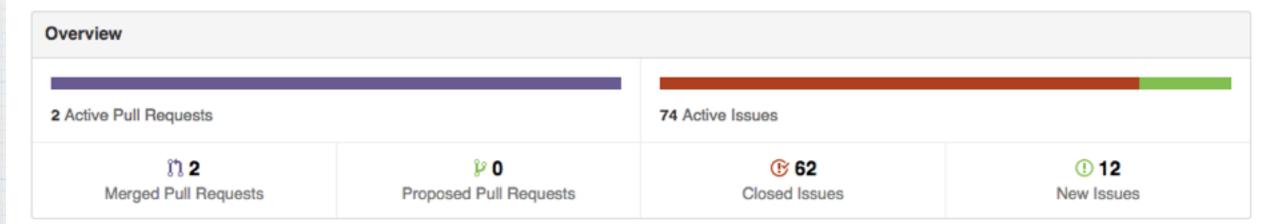

Excluding merges, 5 authors have pushed 66 commits to master and 66 commits to all branches. On master, 68 files have changed and there have been 1,085 additions and 393 deletions.



small RNA-seq analysis

processing & QC

cutadapt fastqc qualimap multiqc detection & annotation

miraligner tdrmapper

de-novo

seqcluster mirdeep2 for mirna protac for pirna (next)

detection & annotation

trimmed and collapsed reads


miraligner (miRBase)

tdrMapper (tRNA reference)




isomiRs package

isomiRs at 5' end of the miRNAs

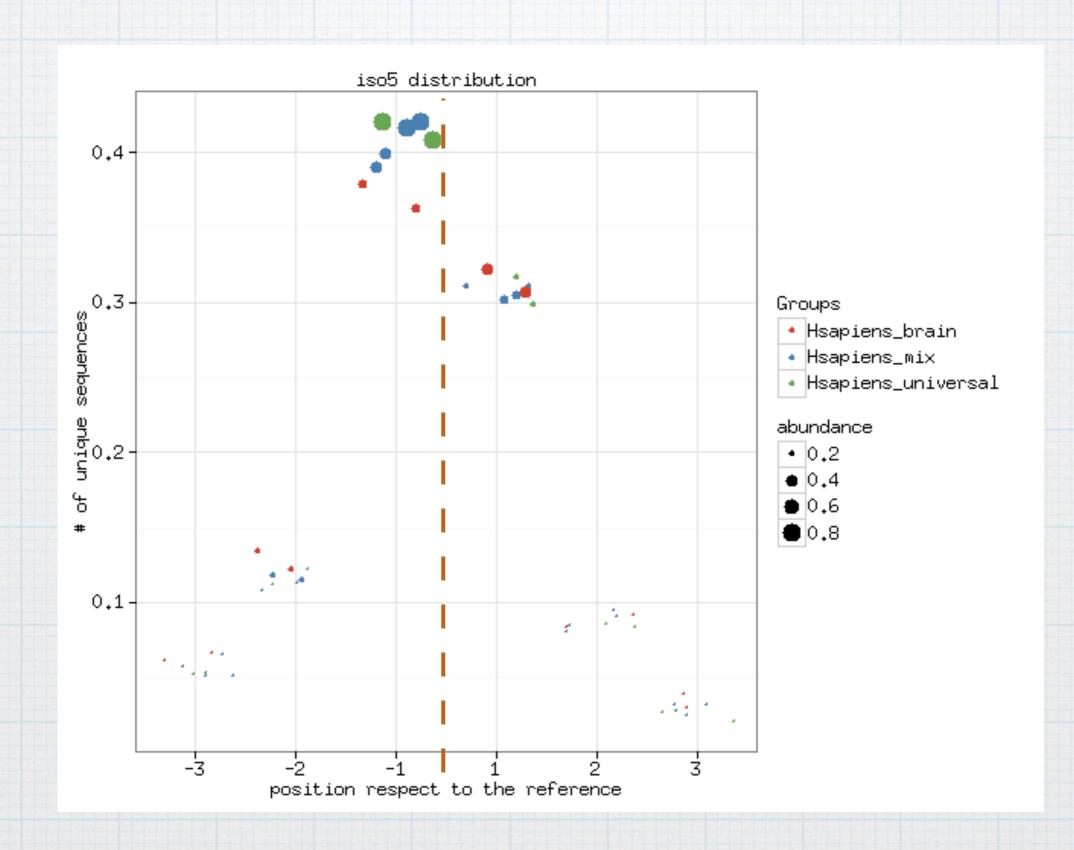

tRNA analysis

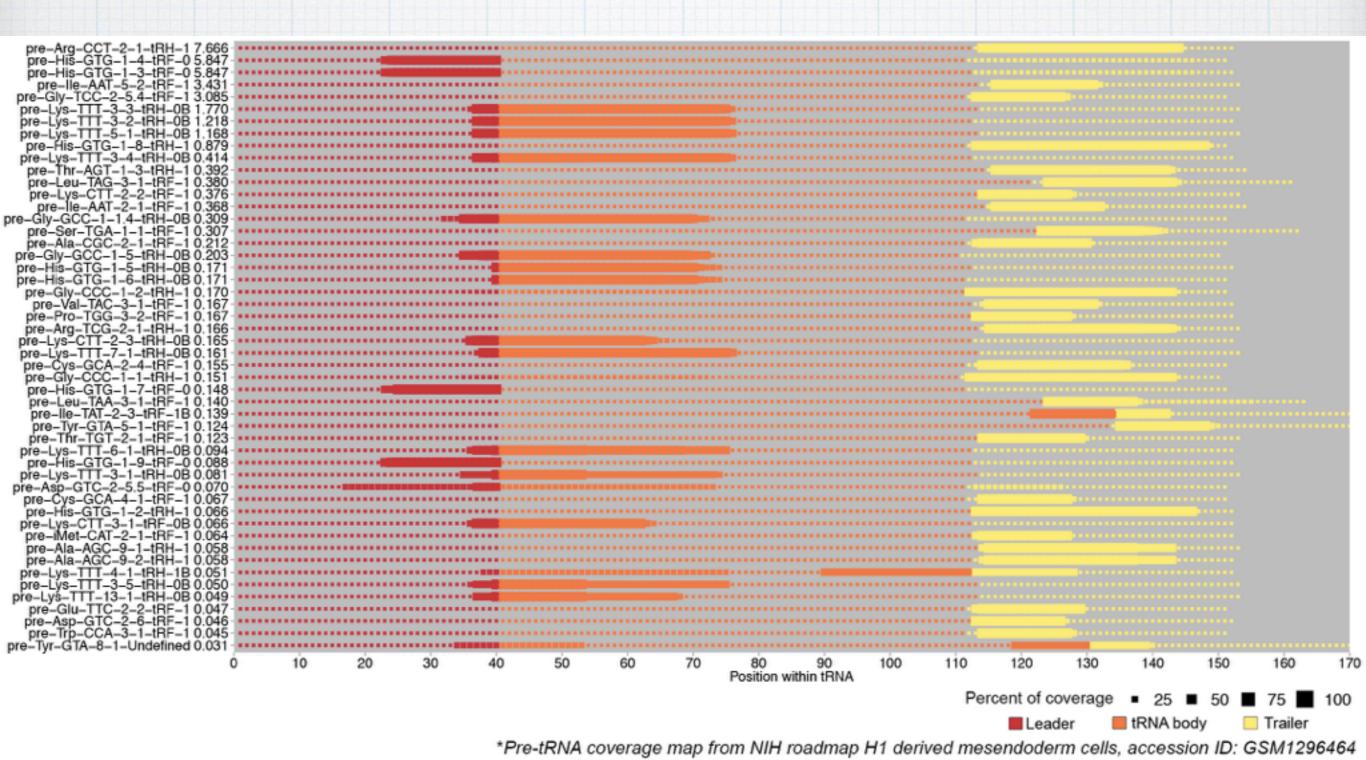

de-novo detection

trimmed and collapsed reads

collapsing samples into one

seqcluster (genome and annotation)


align to genome


mirdeep2 (genome and annotation)


visualization & quantification

expression matrix

miRPeep2 output

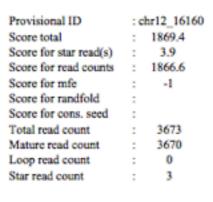

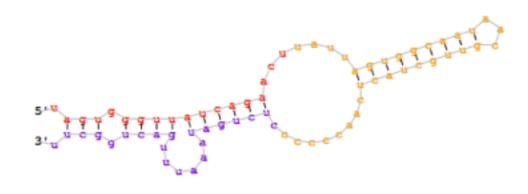

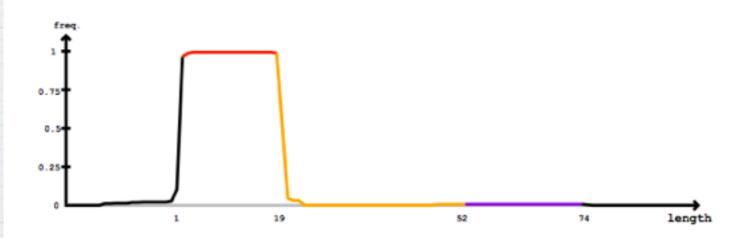

	Mature	Star			
5'-	gaaaugocugguocaaugguaguggguuaucagaacu:auuaguggcaauaaoguugcuacucaacoc	ccucugauanauuugacuggcuuaaaaaanaaaaaaaaucaggaa	-3' obs		
	gaaaugocugguccaaugguaguggguuaucagaacusauuaguggcaauaacguugcuacucaacoc	cocucuga ua a a uuuga cuggo uu aa aa aa aa aa aa aa aa acaggaa	exp		
	(((((()))))))))))))))))	reads	mm	sample
	gocugguccGaugguaguggguuaucagaacu		38	1	seq
	gocugguccGaugguaguggguuaucagaacuua		3	1	seq
	gocugguccGaugguaguggguuaucagaacuuauu		3	1	peq
	cugguccGaugguaguggguuaucagaacuu		3	1	seq
	cugguccGaugguaguggguuaucagaacuua		3	1	seq
	cugguccGaugguaguggguuaucagaacuuauu		3	1	seq
	guccGaugguaguggguuaucagaacuu		16	1	seq


segcluster




meta-cluster

seqcluster deals with multi-mapped reads

Step 1: clustering

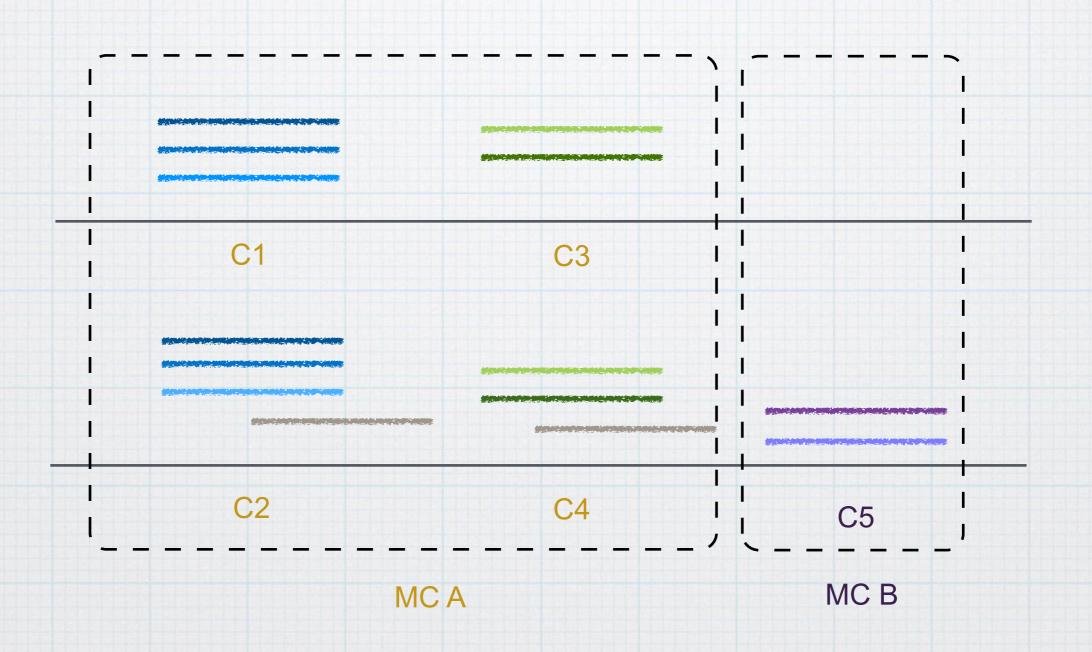

Step 2: cleaning

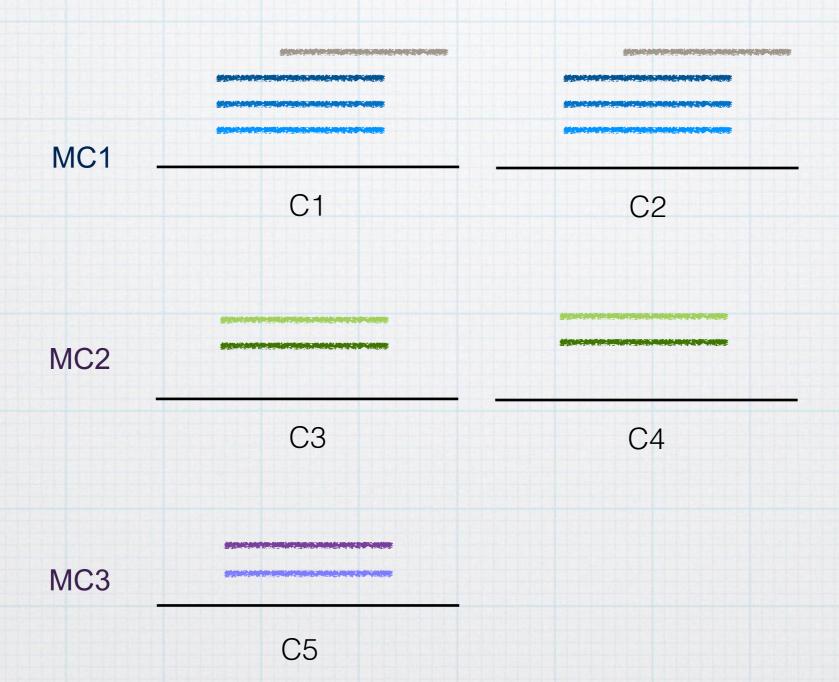

seacluster visualization

	clusters information +						₹
<	> C iii file:///Users/lpantano/repos/seqcluster	rViz/reader.html				Ψ	0 0 0 0 <u>0</u>
+	🗾 BioC 3.3: BUILD/CHE 🏃 🗥 Files - OneDrive 💝 Home - Dr	ropbox H Timesheet - HBC - Ha	Niner Cogalicious - M	 Build your own images 	S Home - PubMed - NO	Google Scholar	
+	Browse						
	Clusters Filter:	Clusters Id:					
	Table with clusters			Table with Locus			
	Sel. I.D. D	Description:		I.D.	Index	Locus:	
	Abundance profile along precursor			Secondary structu	re		

https://github.com/lpantano/seqclusterViz

MultiQ




Phil Ewels ewels


Bioinformatician working with next generation sequencing data.

- Science for Life Laboratory
- Stockholm, Sweden
- phil.ewels@scilifelab.se
- http://phil.ewels.co.uk
- (b) Joined on Nov 3, 2010

48 Followers 21 Starred 23 Following


STAR: % Uniquely mapped reads

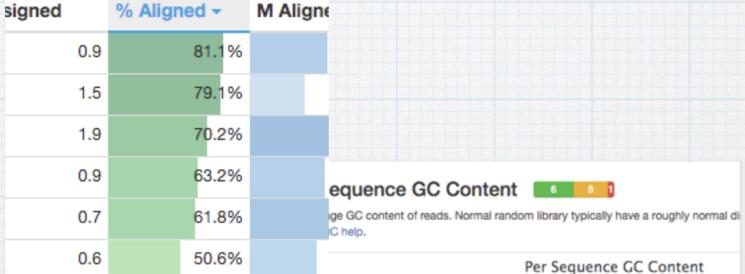

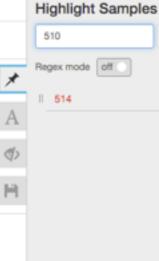

MultiQC Toolbox

.

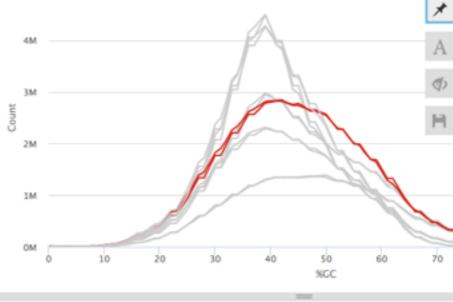


mikQC project

Evaluation of quantitative miRNA expression platforms in the microRNA quality control (miRQC) study

Affiliations | Contributions | Corresponding author

Nature Methods 11, 809–815 (2014) | doi:10.1038/nmeth.3014
Received 27 February 2014 | Accepted 22 May 2014 | Published online 29 June 2014 |
Corrected online 30 July 2014


Quality Control samples

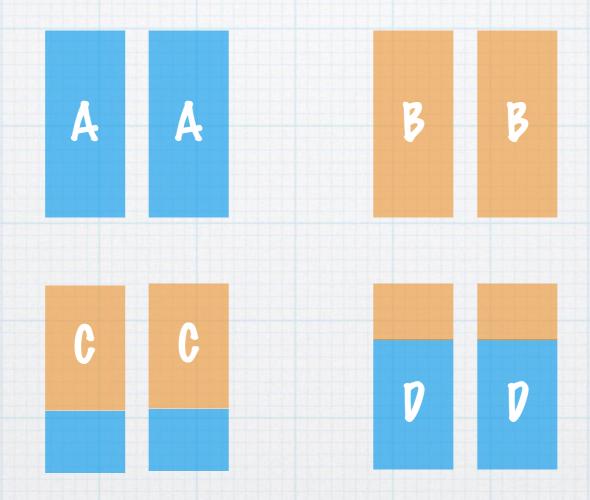

For each molecule:


- * If A > B then A > D > C > B
- * If B > A then A < D < C < B

Good samples

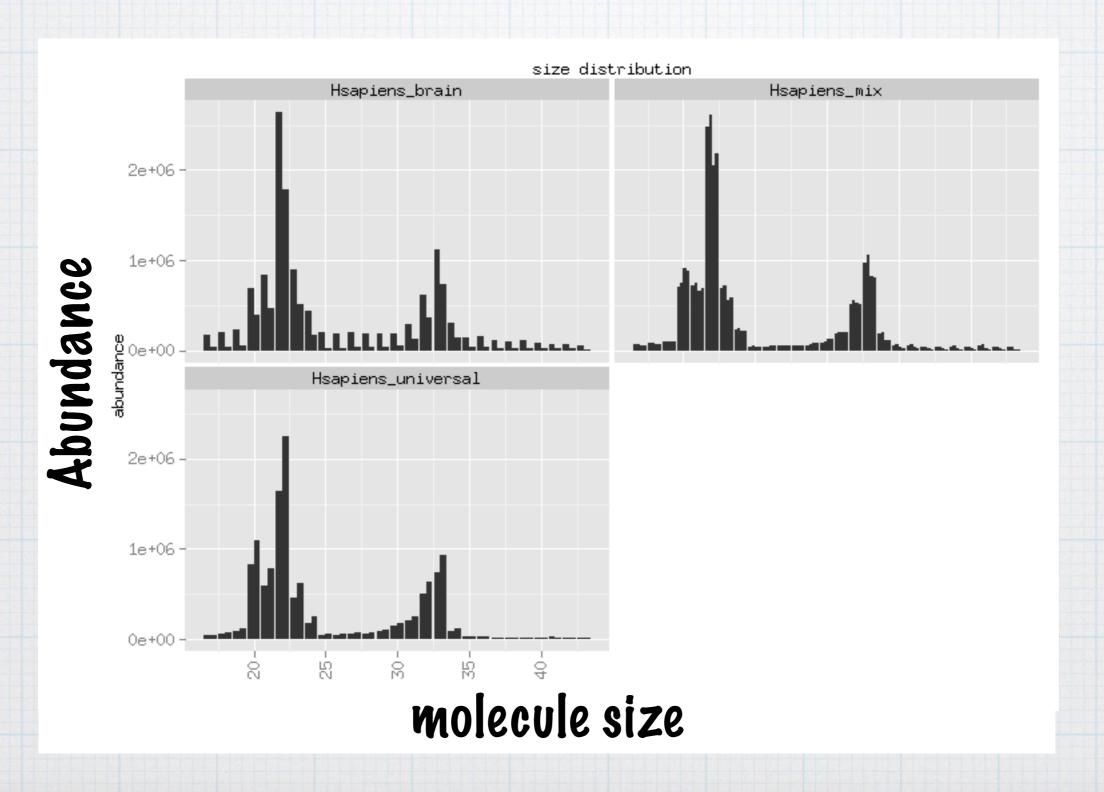


miRNA quantification

miRNAs > 5 counts in average upper quantile normalization

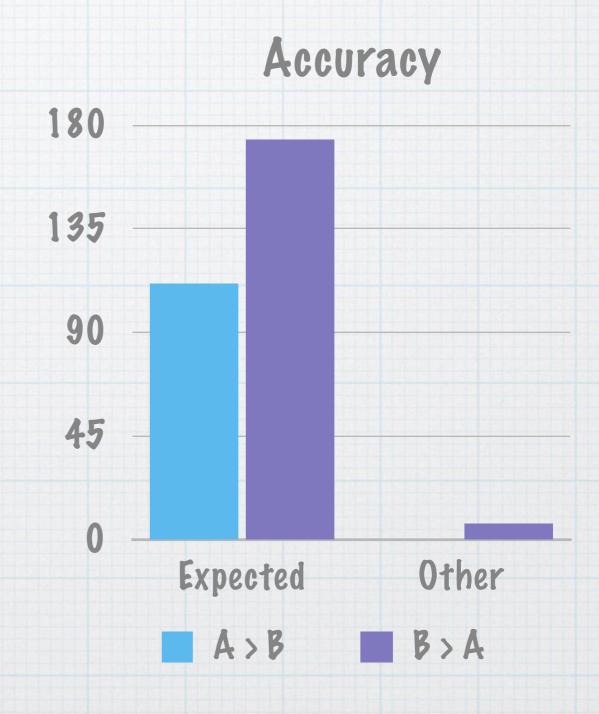


clusters quantification

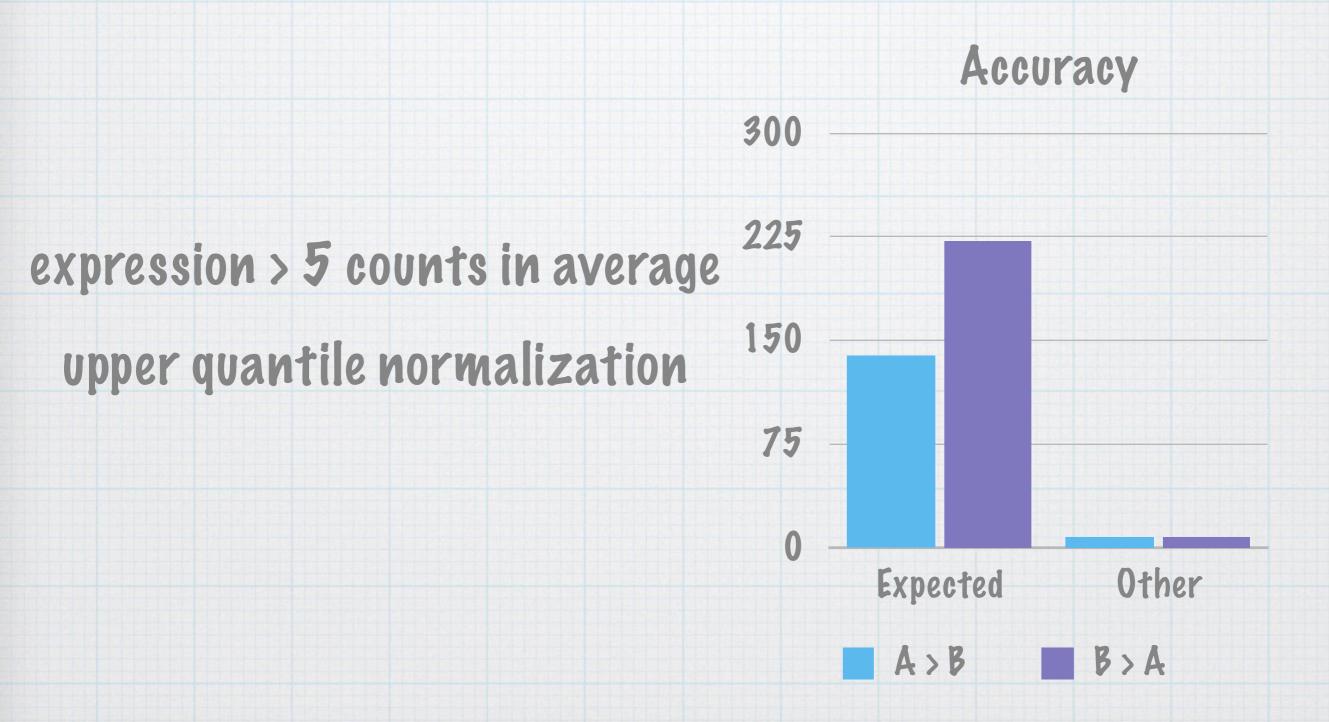

Positive controls

Specificity

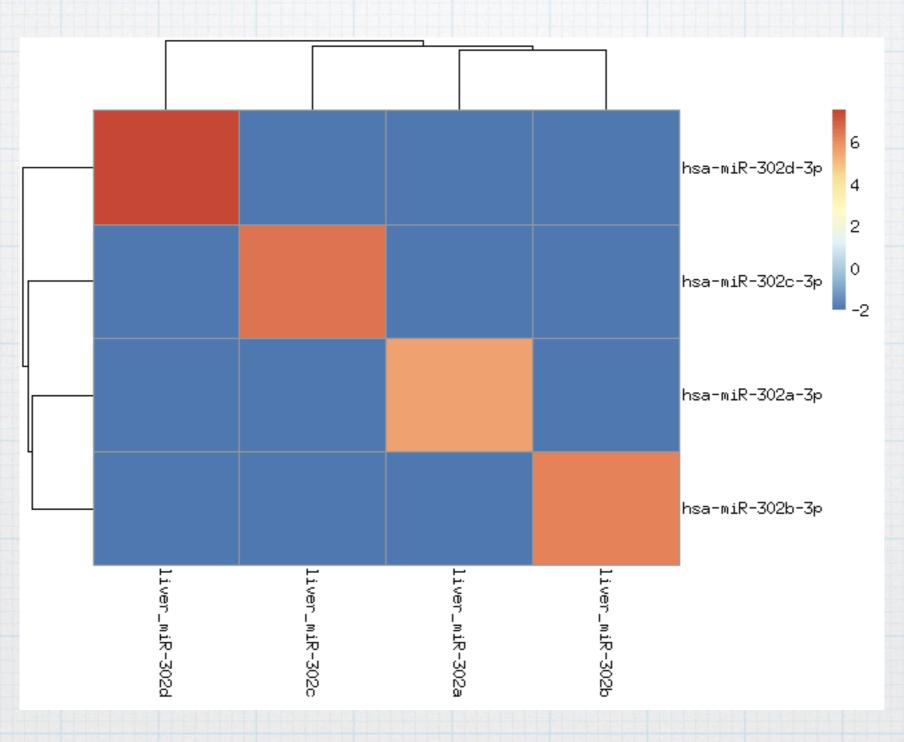

Resources

	Time (h)
organize	0:01
adapter	0:27
alignment	0:26
annotation	3:43
cluster + mirdeep2	4:15
qc	0:04

The time for 8 samples with 6 millions reads each was 8 hours and 57 minutes.

open project for small RNA annotation and analysis



incubator

mirtop

Filters v

Where all ideas and discussions happen to lead to new repositories

Updated 3 days ago

standard formats naming rules

best-practices

New repository

#1 P1

Python # 0 P 0

command lines tool to annotate miRNAs with a standard mirna/isomir naming Updated 3 days ago

miRNAs, tRNAs ...

miRTOP.github.io

project for small RNA standard annotations

Updated on Mar 29

CSS #0 PO

thanks

- * Harvard T.H. Chan School of Public Health
- * Research Computing at Harvard Medical School: Chris Botka, Director of Research Computing and all the people in the team.
- * Special thanks to the authors of those papers to make data available.